BBRYANCE CLEANPIC+

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BBRYANCE CLEANPIC+

User manual

Model : BBPIC001 DENTAL CALCULUS REMOVER CLEANPIC+
Voltage : 3,7V
Power : 1W
Batch number : PIC202004
 Made in China


CONGRATULATIONS ON YOUR PURCHASE !

Before using this product, read carefully this user manual to learn about features and recommendations.
 Keep this notice for future use.

Note: Fully charge the BBRYANCE CleanPic+ toothbrush prior to using it for the first time.

PRODUCT STRUCTURE

With its small size, this product is easy to carry. It's a dental cleansing device that can be used at home. The dental pic head offers the same structure than those used by dentists. It can effectively remove dental calculus, and also has certain cleaning effect on dental plaque or tooth stains.

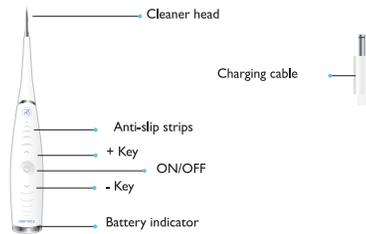PRODUCT FEATURES

- Professional grade alloy metal cleaning head.
- 5 levels of intensity to suit your needs.
- Intelligent high frequency vibration, effective removal of dental calculus and stains.
- Smart and efficient design.
- The battery is durable.

CAUTION

- The metal head of the tooth cleaner is not replaceable.
- Do not place product in direct sunlight. Keep in a -10°C to 40°C environment.
- Do not wash or soak the product for a long time.
- This product is restricted to adult and children aged 12 and more. Store product properly after use. Keep out of reach of children.
- This product contains a battery. Do not place it in a fire or a source of heat to avoid fire or explosion.
- In case of product quality problem, please do not disassemble. Return device to customer service.

PRODUCT SPECIFICATIONS

Rated voltage : 3,7V	Rated current : 130mA
Battery capacity : 300mAh	Standby current : $\leq 20\mu A$
Product size : 183.6x30x22.3mm	Product weight : 53.6g
Charging voltage : 5V	Motor noise : $\leq 60dB$

USAGE METHOD

Before using it for the first time, charge device completely for 2 hours. Always use device until it is completely discharged before charging again (pic 1).
 Press the power button to turn on. Adjust the intensity of vibrations by pressing "+" or "-" . Start cleaning the part of your teeth where you wish to remove calculus and stains (pic 2).
 After use, rinse the metal head with clean water. Wash other metal parts and wipe the body of the device.
 Use: up to once a day, for one minute. In addition to your usual tooth brushing. Users usually use this cleaner 1 to 3 times per week to maintain clean and calculus-free teeth.

CHARGER

1. This product uses a DC interface charging mode. Take the DC interface charging line and plug it into the back of the device. Insert the USB part of the cable into a charger.
2. When charging, the device locks itself and stops working.
3. In charging mode, the light at the bottom is red. The light is green when fully charged and ready for use (pic 3).

SUPPORT AND WARRANTY

BBRYANCE PIC+ is guaranteed for a 2-year duration from the date of purchase.
 Manufacturing or material defects will be repaired or replaced by BBRYANCE in the condition that an official proof of purchase in the qualifying period is provided.
 Contact BBRYANCE customer service on our website: www.bbryance.fr/contact or directly by phone.
WARRANTY EXCLUSIONS
 What is not covered by warranty:
 - Dental pick heads.
 - Damage caused by use of unauthorized replacement parts.
 - Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
 - Normal wear and tear, including chips, scratches, abrasions, or color fading.
 - The proof of purchase has been altered, defaced or removed.
 All implied warranties, including warranties of merchantability and fitness for a particular purpose are limited to the duration of this warranty.
 BBRYANCE or any of its affiliated or subsidiary companies shall not be liable for any special, incidental, consequential or other damages based upon breach of warranty, breach of contract, negligence, tort or any other legal theory.
 This warranty does not give right to claim for compensation.
LEGAL WARRANTY OF CONFORMITY (article of the French consumer code)
 Art. L217-4 "The seller must deliver goods which are in conformity with

the sales contract. The seller is also responsible for defects in conformity that result from packaging, and from set-up or installation instructions when such conformity was covered by the contract or was performed under its responsibility".
 Art. L217-12 "Action resulting from defects in conformity lapses after two years from the date of delivery of the article."
 Goods are in conformity with contract:
 1) If appropriate for the use normally expected for such an article and, when necessary:
 (a) comply with the description given by the seller and possess the qualities of the goods which the seller has held out to the consumer as a sample or model
 (b) have the qualities that a purchaser might legitimately expect in view of public statements made by the seller, the producer or his/her representative in advertising or labelling
 2) have the features defined by mutual agreement between the parties or be suitable for any special requirement of the buyer which was made known to the seller and which the latter agreed to.
 Art. 1641 "A seller is bound to a warranty on account of the latent defects of the thing sold which render it unfit for the use for which it was intended, or that diminish it in a way that the buyer would not have acquired it, or would only have given a lesser price for it, had he known of them".
 Art. 1648 "The action resulting from latent defects must be brought by the purchaser within two years after the discovery of the defect".
 If the BBRYANCE PIC+ device is required to be sent back to the retailer or dealer for repair or replacement, the following information must be provided:
 1) Proof of purchase: a copy or the original of your payment receipt, and the order number for orders made online. Make sure that the date of purchase and product model number are clearly visible. If original documents are provided, they cannot be returned.
 2) Send the product that needs repair or replacement: pack product securely to prevent damage in transit. If possible, use the original packing material.
 3) Enclose a description of the breakdown of malfunction. Include your full name, address where the product must be sent back, phone number, email address (required for shipment). Please fill out the form that customer service will have sent to you following your request, and return it to:
BBRYANCE 143 rue d'Ollioules 83140 Six-Fours les Plages FRANCE.