BBRYANCE **WATERPULSE+**

Manuel d'utilisation / User manual

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Contents

Contents	01
UL Standard	02
Product Description	05
Product Accessories	06
Operation Steps	06-09
Mode descriptions	10
Recommended Techniques	11
Disassembly & Assembly Example	12
Battery disposal and maintenance	13
FAQs and solutions	14
After-sales Service	15

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGERS:

To reduce the risk of electrocution:

- ◆ Do not handle charger with wet hands.
- ◆ Do not immerse in water or other liquid.
- ◆ Do not use while bathing.
- ◆ Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- ◆ Do not reach for a product that has fallen into water.
Unplug immediately.
- ◆ Check the charger cord for damage before the first use and during the life of the product.

To reduce the risk of burns, electric shock, fire, serious injury or damage issues:

Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.

Do not use this product if it has a broken cord or plug. In case of not working properly, or it has been dropped into any liquid, contact your local distributor.

Do not direct water flow under the tongue, into the ear, nose or other delicate organs. This product is capable of producing pressures that may cause serious damage to these organs. Refer to operating instructions for correct usage.

Use this product only as indicated in these instructions or suggestions by your dental professional.

It's normal phenomenon if gum bleed a little when using water flosser for the first time. In case of excessive bleeding, please stop using and ask for your dentist

Fill reservoir with water or other solutions recommended by dental professional.

Do not drop or insert any foreign object into the product.

Keep charger and cord away from heated surfaces.

Do not operate where oxygen or aerosol sprays are being used.

Do not use mouth washes/rinses that contain iodine, bleach, or tea treeoil.

Do not use while wearing any oral jewelry. Please take off prior to use.

Do not use if you have an open wound on your tongue or in your mouth.


If your physician advise you to receive antibiotic premedication, you should consult your dentist before using this instrument or any other oral hygiene aid.

For household use only. Do not use outdoors.


Closely instruct and supervise children and individuals with special needs in the proper use of this product.

Children should be supervised to ensure that they do not play with the appliance.

Product Description

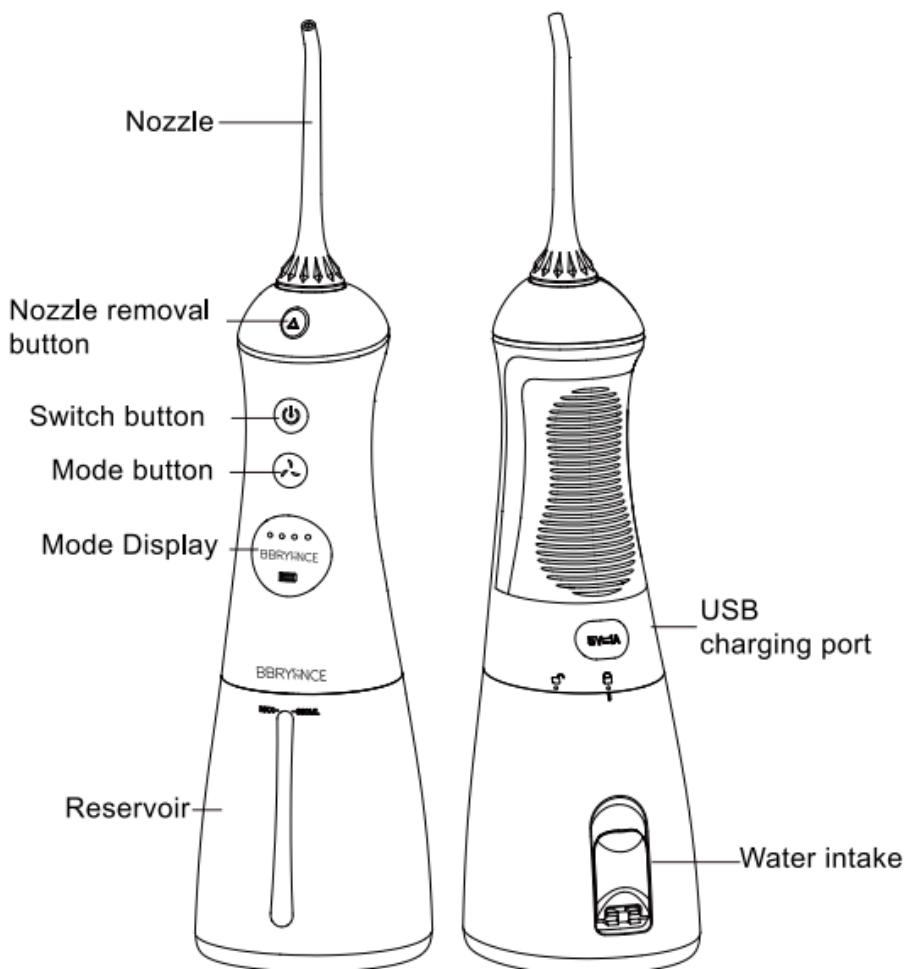Product Accessories

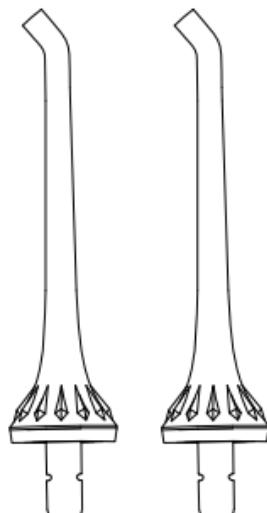Nozzles x 2


Charging cable x1

Replaceable nozzles/accessories can be purchased from the vendor.

Operation steps


Activation

Very important! Please fully charge the device before using it for the first time. Plug the USB cable into the charging socket, and turn on the power. The product cannot be used while charging.

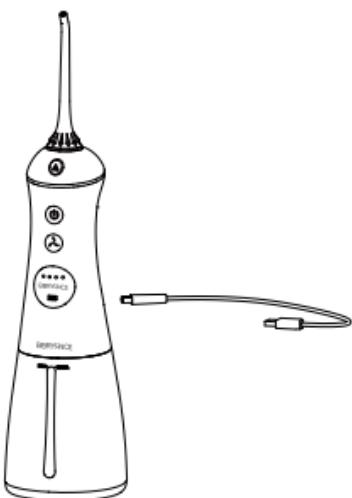
Nozzle installation and removal

Insert the nozzle into the middle hole of the host until it clicks, and the installation is successful. To replace or remove the nozzle, please press and hold the nozzle removal button and then take out the nozzle.

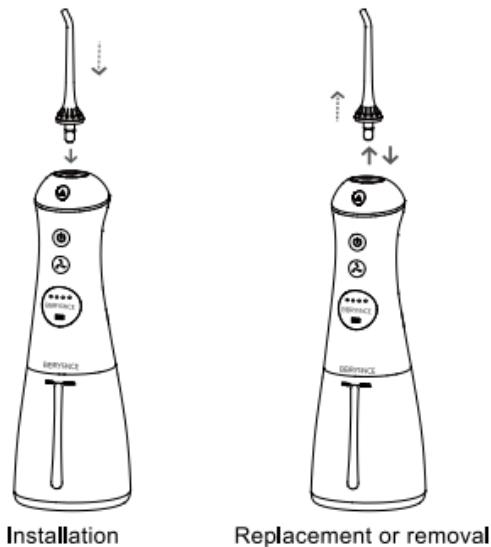As shown in Figure 1, open the reservoir cover and then pour water into it;

Note: Please make sure that the device is turned off when filling the reservoir.

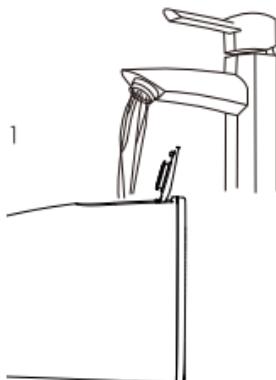
- ◆ The power switch and the mode buttons are on the panel. the power button controls device activation and deactivation, and the mode button is responsible for switching between the 6 modes.
- ◆ If the switch button is pressed again when the product is in a working state, the product will stop working. The regular working time of the product is 2 minutes, and if the product is to be used again, press the power button again.
- ◆ When charging, a indicator lamp will flash in orange. When it stops flashing in white, the product is fully charged. When the battery is low, the light flashes to remind you to recharge. When the battery is too low, the device's function cannot be activated by pressing the power button, and the approximate battery power is displayed when it is turned on; when it is turned off, the approximate battery power of the product can be detected by pressing the mode button (red, blue, and green).

- ◆ The charger is suitable for: AC 100-240VAC, 50/60Hz with an output of 5V 1A.
- ◆ The device cannot be used when charging.
- ◆ When in use, keep it upright to ensure that the equipment is in the best working condition; by rotating the nozzle, the direction of the water flow can be adjusted.
- ◆ Gum bleeding may occur the first time you use the oral irrigator, because everyone can withstand different pressures in their gums. This is normal, so please rest assured!
- ◆ When using the product for the first time, please start at the lowest gear and gradually increase the pressure as your oral cavity adapts.

If the water tank cover falls off

1. Align the gasket of the water tank cover with the water tank mouth;
2. Press it in until you hear a "click" sound.

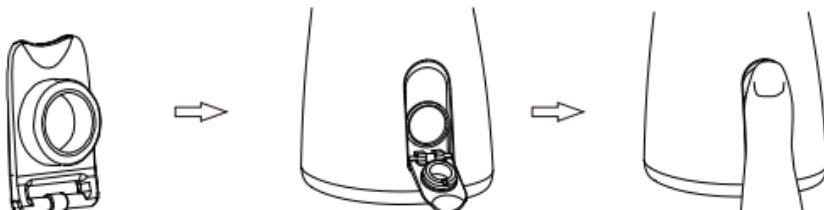Mode descriptions

Gear 1: When 1 white lamp is always on, it signals the 1st gear, and the water pressure is at 10-40PSI.

Gear 2: When 2 white lamps are always on, it signals the 2nd gear, and the water pressure is at 30-60PSI.

Gear 3: When 3 white lamps are always on, it signals the 3rd gear, and the water pressure is at 50-90PSI.

Gear 4: When all 4 white lamps are always on, it signals the 4th gear, and the water pressure is at 60-110PSI.

Gear 5: When the 4 white lamps flash from left to right, it signals the pulse gear, and the water pressure is at 60-110PSI.

Gear 6: When the LOGO is flashing with a pulsating light, it signals the DIY gear, and the water pressure is at 10-110PSI.


Setting the DIY gear: In any mode, press and hold the mode button for 1 second, and the logo will flash rapidly; quickly press the mode button to increase the water pressure by one level (for DIY gears, 8 levels of water pressure can be set, where every flashing lamp signals 1 level of water pressure, and every lamp that is always on signals two levels of water pressure). Press and hold the mode button for 1 second or do nothing for 5 seconds to exit the setting mode.



Recommended Techniques

Lean over the wash basin and put the tip in the mouth. Aim the tip at your teeth.

Turn on the switch and shift to the suitable mode.

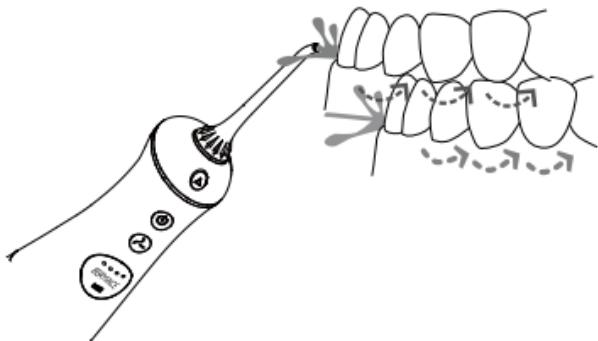
Make water jet at 90° to the gingiva, and close your lips slightly to avoid water splash everywhere. For best effect, start with the molars and clean towards the front teeth. Slide the tip along the gumline and keep flossing a while at inter-dental place, till everywhere is cleaned.

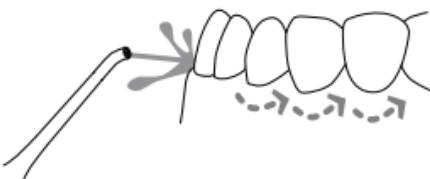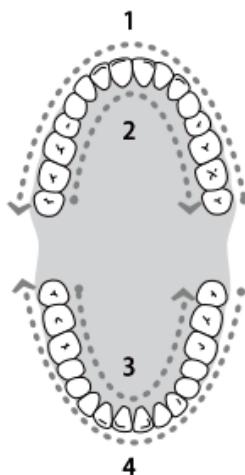
Disassembly & Assembly Example

When Finished

Turn the unit off. Lift the flip top on reservoir as shown, and empty any liquid left in the reservoir. Turn the unit on till reservoir is empty and turn the unit off.

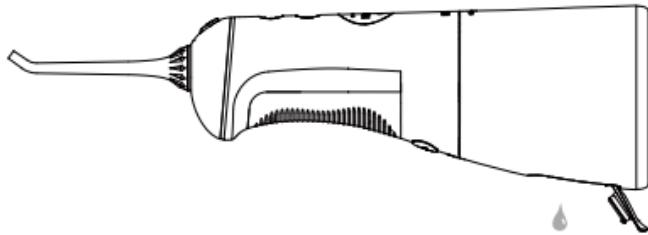To remove the reservoir, hold the unit with your left hand and slide the reservoir down counterclockwise with your right hand.

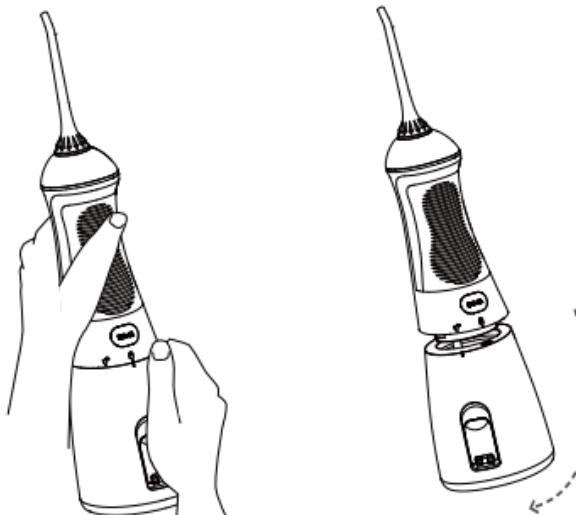Note: There may be moisture inside the tank (tank and suction pipe), which is the residual distilled water used for product testing.

Battery Disposal and maintenance

Battery maintenance


- ◆ The device must be fully charged before first using, and subsequent charges take about 2 hours;
- ◆ Plug the USB cable into the charging socket and connect it to a power source;
- ◆ If the device will be stored for a long time, please fully charge the device before storing;
- ◆ In order to maximize the service life of the battery, please charge it immediately when the battery is low;

Battery handling


This device contains a lithium battery, which is to be recycled at the end of its life.

Please recycle the battery at a local recycling center. The battery cannot be replaced.

Removal of the built-in battery

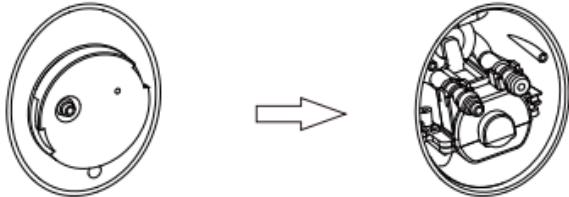Remove the nozzle and water tank, and pry open the base with a screwdriver or any other tool.

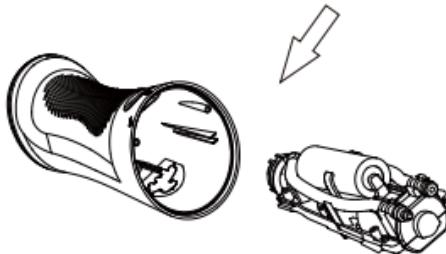Remove the battery for recycling

Be careful not to short-circuit the positive and negative ends of the removed battery, and insulate the ends with tape.

FAQs and solutions

Problems	Reasons	Solutions
Leakage between the nozzle and the body of the device	The seal ring in the nozzle holder is damaged, the nozzle is not installed securely	Contact the vendor for repair, install the nozzle securely
Leakage in the water pump inlet	The seal ring is damaged, or there are impurities on the inlet side	Contact the vendor for repair
Leaking water tank	Damaged seal	Replace the seal ring, contact the vendor for repair
Not enough pressure	Low voltage in use	Use after fully charged
Device cannot start	Low voltage in use	Use after fully charged or contact the vendor for repair
The product does not spray water	First use or not used for a long time	After adding a proper amount of water to the nozzle installation port, reinstall the nozzle, and then turn on the power switch

If you still can't solve the problem according to the above methods, please consult the vendor shop or go to a designated after-sales repair point for repair.

After-sales service

BBRYANCE WATERPULSE water flosser is guaranteed for a 24-month duration from the date of purchase. Manufacturing or material defects will be repaired or replaced by BBRYANCE in the condition that an official proof of purchase in the qualifying period is provided.

For warranty service, contact BBRYANCE customer service on our website: www.bbryance.fr/contact or directly by phone.

WARRANTY EXCLUSIONS

What is not covered by warranty:

- Nozzles, charging cable.
- Damage caused by use of unauthorized replacement parts or unsuitable brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, or color fading.
- The proof of purchase has been altered, defaced or removed.

All implied warranties, including warranties of merchantability and fitness for a particular purpose are limited to the duration of this warranty. BBRYANCE or any of its affiliated or subsidiary companies shall not be liable for any special, incidental, consequential or other damages based upon breach of warranty, breach of contract, negligence, tort or any other legal theory.

This warranty does not give right to claim for compensation.

LEGAL WARRANTY OF CONFORMITY (article of the French consumer code)

Art. L217-4 "The seller must deliver goods which are in conformity with the sales contract. The seller is also responsible for defects in conformity that result from packaging, and from set-up or installation instructions when such conformity was covered by the contract or was performed under its responsibility".

Art. L217-12 "Action resulting from defects in conformity lapses after two years from the date of delivery of the article."

Goods are in conformity with contract:

1) If appropriate for the use normally expected for such an article and, when necessary:

(a) comply with the description given by the seller and possess the qualities of the goods which the seller has held out to the consumer as a sample or model

(b) have the qualities that a purchaser might legitimately expect in view of public statements made by the seller, the producer or his/her representative in advertising or labelling

2) have the features defined by mutual agreement between the parties or be suitable for any special requirement of the buyer which was made known to the seller and which the latter agreed to.

Art. L641 "A seller is bound to a warranty on account of the latent defects of the thing sold which render it unfit for the use for which it was intended, or that diminish it in a way that the buyer would not have acquired it, or would only have given a lesser price for it, had he known of them".

Art. L648 "The action resulting from latent defects must be brought by the purchaser within two years after the discovery of the defect".

If the BBRYANCEWATERPULSE device is required to be sent back to the retailer or dealer for repair or replacement, the following information must be provided:

1) Proof of purchase: a copy or the original of your payment receipt, and the order number for orders made online. Make sure that the date of purchase and product model number are clearly visible. If original documents are provided, they cannot be returned.

2) Send the product that needs repair or replacement: pack product securely to prevent damage in transit. If possible, use the original packing material.

3) Enclose a description of the breakdown of malfunction. **Include your full name, address where the product must be sent back, phone number, email address** (required for shipment). Please fill out the form that customer service will have sent to you following your request, and return it to:

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